

Reduce Spinout in Viscous Solutions!

Rare Earth – Samarium Cobalt Spinbar® Magnetic Stirring Bars


Reduce frequency of spinout in viscous solutions or when stirring at high speeds. The high magnetic field strength of Samarium Cobalt (Rare Earth) magnets can provide up to three times stronger coupling vs. Alnico magnetic stirring bars.

- Stronger coupling with the drive magnet keeps the stirring bar centered and rotating effectively
- Encapsulated in bright green, FDA grade, USP Class VI Teflon® PTFE; easy to differentiate from other stirring bars
- Autoclave safe
- Manufactured in a registered ISO 9001:2015 facility, as verified by SGS Certification

**3X Stronger Stirring Power*
with Samarium Cobalt**



Samarium Cobalt vs. Alnico Magnetic Stirring Bar
Revolutions per Minute without Spinout in Fluids of Increasing Viscosity

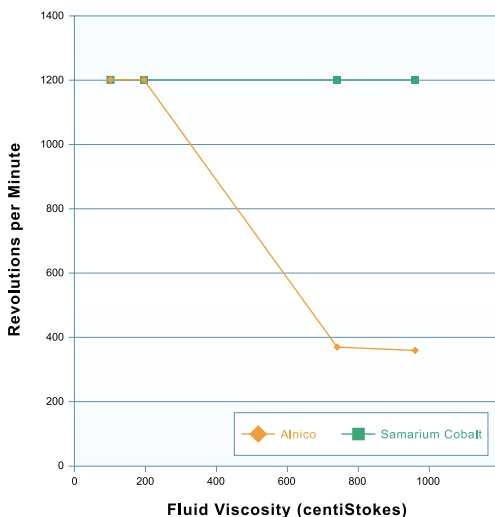

*Testing done with 1/2" x 3/8"H magnetic stirring bars

CATALOG NO.	SHAPE	DIMENSIONS (L X DIA.)
F37102-1238	Octagon	1/2 x 3/8" (1.30 x 0.95cm)
F37102-0138	Octagon	1 x 3/8" (2.54 x 0.95cm)
F37102-0238	Octagon	2 x 3/8" (5.10 x 0.95cm)
F37102-0003	Octagon	3 x 1/2" (7.60 x 1.30cm)
F37174-0000	Fluted	2 x 8/10" (50 x 21mm)

Teflon® is a registered trademark of The Chemours Company FC, LLC

www.belart.com • 1-800-423-5278

661 Route 23 South, Wayne, NJ 07470