


INSTRUCTIONS



Drying Tubes With Tube Fittings

Catalog Number: 19960, 61, 62-0000
Extra Fittings Cat. No. 19965-0000

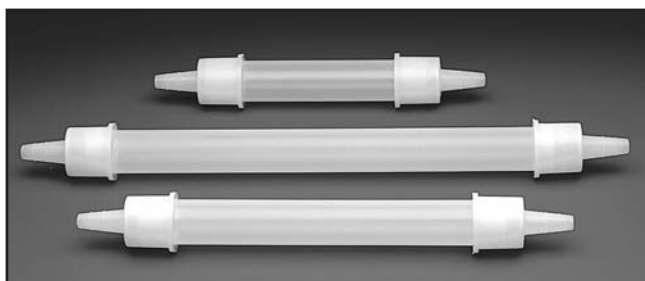

The drying tubes are made of polyethylene and should not be placed in an autoclave. The end fittings are molded of polypropylene and are autoclavable at 250°F (121°C). The tubes are supplied empty and may be filled with a drying agent. The drying tubes may vary slightly in diameter. For this reason, the end fittings are

molded with a "stepped" inner diameter. If the tubes are on the upper limit of the diameter ranges, the fitting will slip on only about 1/4" (6.4mm). It is not necessary for the tube to enter the fitting completely. If the tubes are on the lower limit of the diameter range, they will fully enter the fitting.

BEL-ART PRODUCTS

661 Route 23 South, Wayne, NJ 07470-6814 USA
TEL: 1-800-4BEL-ART • FAX: 973-694-7199 • www.belart.com

919650002

Scienceware® is a registered TM of Bel-Art Products
Bel-Art Products assumes no obligation or liability for any advice furnished by it, or for results obtained with respect to these products. All such advice is given and accepted at the buyer's risk.
©2008 Bel-Art Products. All Rights Reserved.


03/12

INSTRUCTIONS



Drying Tubes With Tube Fittings

Catalog Number: 19960, 61, 62-0000
Extra Fittings Cat. No. 19965-0000

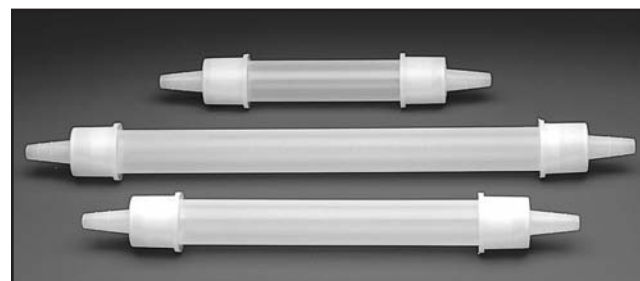

The drying tubes are made of polyethylene and should not be placed in an autoclave. The end fittings are molded of polypropylene and are autoclavable at 250°F (121°C). The tubes are supplied empty and may be filled with a drying agent. The drying tubes may vary slightly in diameter. For this reason, the end fittings are

molded with a "stepped" inner diameter. If the tubes are on the upper limit of the diameter ranges, the fitting will slip on only about 1/4" (6.4mm). It is not necessary for the tube to enter the fitting completely. If the tubes are on the lower limit of the diameter range, they will fully enter the fitting.

BEL-ART PRODUCTS

661 Route 23 South, Wayne, NJ 07470-6814 USA
TEL: 1-800-4BEL-ART • FAX: 973-694-7199 • www.belart.com

919650002

Scienceware® is a registered TM of Bel-Art Products
Bel-Art Products assumes no obligation or liability for any advice furnished by it, or for results obtained with respect to these products. All such advice is given and accepted at the buyer's risk.
©2008 Bel-Art Products. All Rights Reserved.

03/12